NOOA

PRO PERFORMANCE

# Electronic Personal Scale

NOPS70MB

User Manual

GB

Thank you for purchasing your new NOOA electronic personal scale. These operating instructions will help you use it properly and safely. We recommend that you spend some time reading this instruction manual in order that you fully understand all the operational features it offers. Read all the safety instructions carefully before use and keep this instruction manual for future reference.

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## SAFETY WARNING


1. Read all the instructions carefully before using this product and keep them for future reference.
2. This product has delicate electronic components. Interference may occur when it is near devices which transmit radio signals. If errors occur in the display, remove the product away from these signals or remove the battery briefly and then insert it again.
3. Keep batteries away from children. Do not leave batteries lying around. They could be swallowed by children or pets. In the event of swallowing a battery, contact a doctor immediately.
4. This scale is for household use only. Do not use for commercial purpose.
5. The glass platform will be slippery when it is wet. Keep it dry!
6. Stand still while being weighed.
7. Do not strike, shock or drop the scale.
8. Treat the scale with care as it is a precise instrument.
9. Clean the scale with a damp cloth and avoid water getting inside the scale. Do not use chemical/abrasive cleaning agents.

10. Keep the scale in cool place with dry air.
11. Always keep your scale in horizontal position.
12. If the scales fails to turn on, check to see if the battery is installed, or if the battery power is low. Insert or replace new battery.
13. If there is an error on the display or you can not turn off the scale for a long time, take out the battery for about 3 seconds, then re-install it to remove the soft-fault. If you can not resolve the problem, contact your dealer for consultation.

## PRODUCT OVERVIEW


## TO REPLACE BATTERY

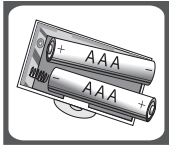
1. Powered by 2x1.5V AAA batteries(not included).
2. Take out used batteries, with the aid of a sharp object if needed. Install new batteries by putting one side of the battery down first and then pressing down the other side.  
(Observe the polarity)

## LCD DISPLAY
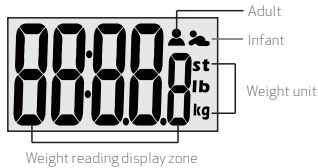
## BUTTON FUNCTION

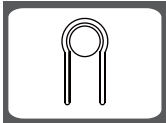Button on scale bottom


Under weighing mode:

1. Long press: Convert mode between Adult mode (normal weighing) and infant mode
2. Short press: Convert unit among kg/lb/st.


Remark: The unit system of the scale depends on the sales destination or the requirement of distributor.

## OPERATION

### 1. Adult Mode - Normal weighing

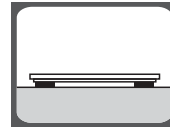
Step 1: To ensure accuracy, place the scale on a hard, flat floor and avoid carpet or soft surface.

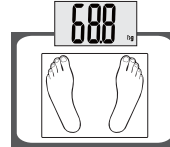Step 2: Gently step on the scale, then the scale will be automatically turned on. Stand evenly on the scale without moving and wait until the LCD shows the weight value. The weight reading will flash three times and locked (as shown on picture).


### 2. Infant Mode - Weighing infant while being held by an adult

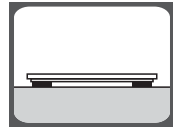Step 1: To ensure accuracy, place the scale on a hard, flat floor and avoid carpet or soft surface.


Step 2: Adult gently steps on the scale, then the scale will automatically turn on. Stand evenly on the scale without moving and wait until the LCD shows the weight value of the adult. The weight reading will flash three times and return to zero. Infant icon will then flash. (as shown on picture)

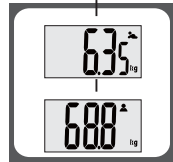
Step 3: Adult steps off the scale and then steps on again holding the infant. Stand evenly on the scale without moving and wait until the LCD shows the weight value of the infant. The weight reading will flash and lock for about 8 seconds. LCD will then change to indicate the weight value of the adult again.

## AUTOMATIC SWITCH OFF

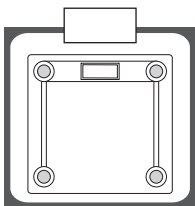When you step off the scale, it will automatically switch off. Auto shut-off occurs if the display shows the same weight reading for approximate 20 seconds.

## WARNING INDICATIONS
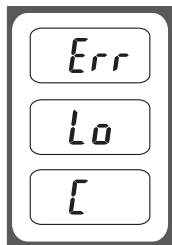
### OVERLOAD INDICATION

Scale overload. (max. weight capacity is indicated on the scale). Step off to avoid damage.

### LOW BATTERY INDICATION

Battery power is low. Install new battery.

### RE-TEST

Display "L" when you step off the scale, means this weight reading has an error.

## Features/Specifications

1. Model: NOPS70MB
2. 74 x 40 mm LCD
3. Low battery / over load indication
4. Capacity: 180kg / 400lb / 28st
5. Graduation: 100g / 0.2 lb (adult mode); 50g / 0.1lb (infant mode)
6. Power: 2 x 1.5V AAA

This symbol on the product or in the instructions means that your electrical and electronic equipment should be disposed at the end of its life separately from your household waste. There are separate collection systems for recycling in the EU.

### Notes for battery disposal

The batteries used with this product contain chemicals that are harmful to the environment. To preserve our environment, dispose of used batteries according to your local laws or regulations. Do not dispose of batteries with normal household waste.

For more information, please contact the local authority or your retailer where you purchased the product.

Imported and Exclusively marketed by:

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