The logo for 'point' is displayed in white lowercase letters on a red rectangular background with rounded corners. A registered trademark symbol (®) is located at the top right of the word.

Renew Body Composition Scale

POPS3000
User Manual

The logo consists of the words 'LIFE', 'MADE', and 'EASIER' stacked vertically. 'LIFE' is in a large, bold, grey font. 'MADE' is in a smaller, bold, orange font. 'EASIER' is in a large, bold, grey font. The entire logo is enclosed in a thin red rectangular border.

GB

Thank you for purchasing your new POINT RENEW body composition scale. These operating instructions will help you use it properly and safely. We recommend that you spend some time reading this instruction manual to fully understand all the operational features it offers. Read all the safety instructions carefully before use and keep this instruction manual for future reference.

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SE

Tack för att du köpt en ny POINT RENEW kroppsanalysvåg. These operating instructions will help you use it properly and safely. We recommend that you spend some time reading this instruction manual in order that you fully understand all the operational features it offers. Read all the safety instructions carefully before use and keep this instruction manual for future reference.

Instruction Manual - English	page	4	-	15
Bruksanvisning- Norsk	sida	16	-	27
Käyttöohjeet - Suomi	side	28	-	39
Brugsanvisning - Dansk	side	40	-	51
Bruksanvisning - Svenska	sivu	52	-	63


SAFETY WARNING

1. Always use the scale on a solid and flat surface. The device is designed to allow auto step-on. Ignore the first reading and only start recording measurements from the second weighing.
2. Remove your shoes and socks, and ensure your feet are clean before using the scale.
3. Step onto the weighing platform gently.
4. For consistency, take your measurements at the same time each day.
5. Measurement results may be misleading after intensive exercise, excessive dieting, or extreme dehydration.
6. Keep the scale away from strong electro-magnetic fields when in use.
7. Children under 10 years old or adults over 100 years old should only use the device in normal weighing mode.
8. The measured data may have some deviations for some individuals, such as:
 - Adults over 70 years old
 - Body building athletes or other occupational athletes.
9. People with the following conditions should not use the device: pregnant women, those with symptoms of edema, or on dialysis treatment.

10. The Fat% and Hydration% content evaluated by the Body Composition Scale is only for reference (not for medical purposes). If your Fat% or Hydration% content is above or below normal levels, please consult your doctor for more advice.
11. If the scale fails to turn on, check whether the battery is installed, or if the battery power is low. Insert or replace with a new battery.
12. Do not use for commercial purposes.
13. The platform can be slippery if you weigh yourself on a wet surface. If the weighing platform surface is glazed, never step on the edge of the scale and always keep your balance during your weighing.
14. People with heart pacemakers or other implanted medical devices should not use the body fat testing function.
15. Keep batteries away from children. Do not leave batteries lying around. They could be swallowed by children or pets. In the event of swallowing a battery, contact a doctor immediately.


PRODUCT OVERVIEW

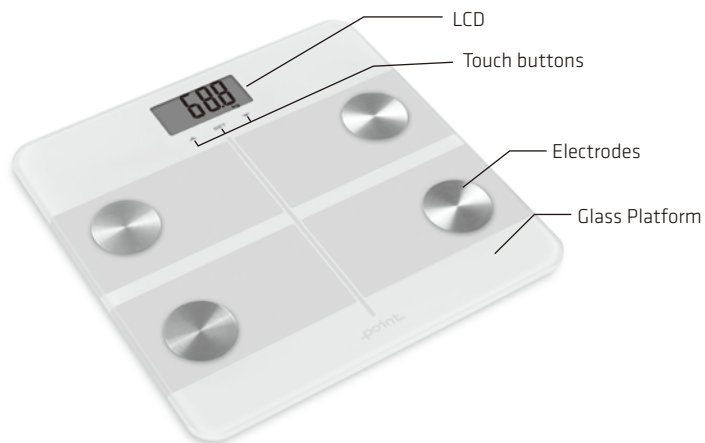BATTERY INSTALLATION


- Take out the used battery with the aid of a pointed object.
- Install a new battery by putting one side of the battery beneath the battery baffle and then pressing down the other side. Please take note of the polarity.

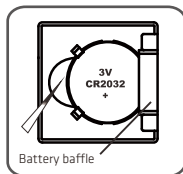KEY ILLUSTRATION


- SET** : Turn on the scale and proceed to the next parameter. (Please refer to the Setting Personal Parameters overleaf)
- ▲ : Increase value with each press.
Press and hold to continuously increase the value.
- ▼ : Decrease value with each press.
Press and hold to continuously decrease the value.

LCD DISPLAY

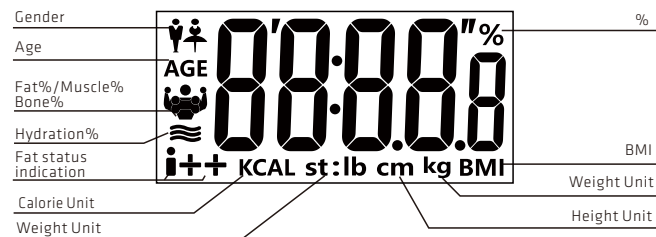TO CONVERT WEIGHT UNIT


- Find the unit switch on the back of the scale.
- Press the unit conversion button when the scale is powered on, and the LCD shows current weight unit. Press the button again to convert the weight unit.

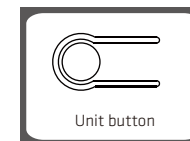NORMAL WEIGHING MODE (by passing parameter setting)

1. Auto step on weighing


After installing the battery, place the scale on a flat floor. Stand still on the scale when testing. Then wait for a few seconds until the display shows the same weight reading for approximately a few seconds.

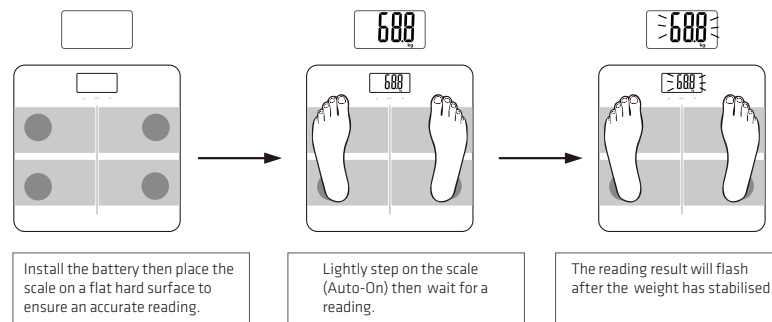2. Consecutive add-on weighing

After the digits have stabilized, you can perform additional weighings. Ensure each add-on loading exceeds 2kg to get a new reading.

For example: If you want to get a weight value of your baby, follow these steps: Weigh yourself first and note the reading. Without stepping off the scale, hold your baby and weigh yourself again. The scale will display the combined weight. Subtract your initial weight from this combined reading to find your baby's weight.

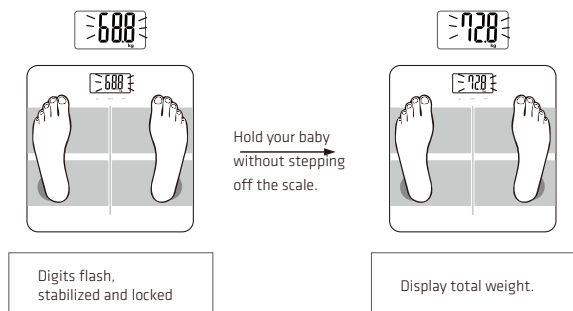FAT/HYDRATION/MUSCLE/BONE/CALORIE MEASURING MODE

(1) Setting personal parameters

A. The first time you install the battery and switch the scale on, all default parameters will be set as follows:


Memory	Gender	Height	Age
1	Female	165cm/5'05.0"	25

B. Parameters can be set in the following range:

Memory	Gender	Height	Age
1~12	male/female	100~250cm 3'03.5"~8'02.5"	10~100

C. Parameter setting

Example: the 3rd user, female, 168cm, 28 years old, can set her parameters as follow:

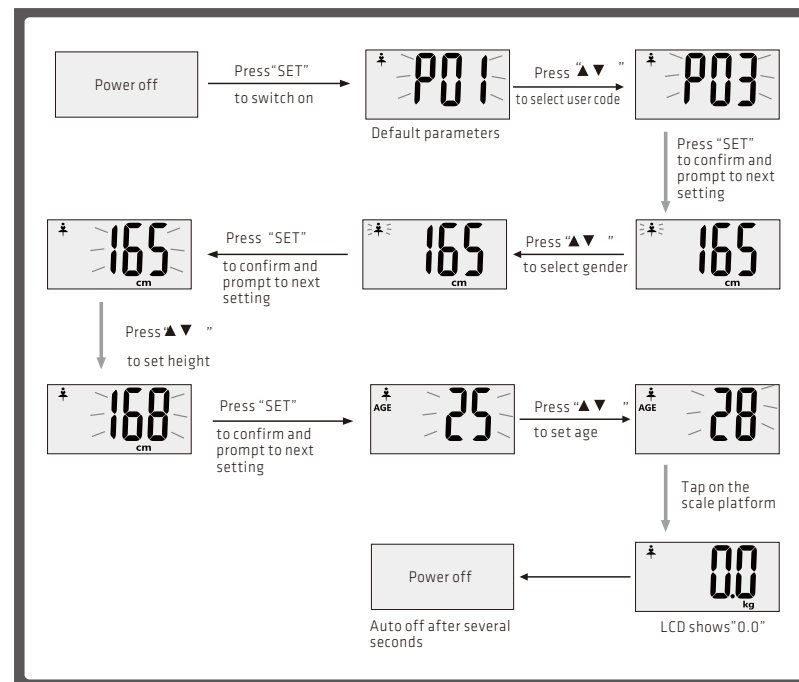
NOTE:

- During parameter setting, if the scale is idle for 6 seconds, the parameters already set will be memorized automatically and the LCD shows "0.0", meaning that body fat measuring will be started.


If the LCD shows "0.0" during the parameter setting, press SET for 3 seconds to return to the setting mode, then you can set the user code, gender, height and age

(2) To Start Measuring (Make sure your feet are touching the electrodes on the scale platform, or your body fat can not be measured.)

To ensure accuracy, please remove your shoes & socks before stepping on.
 Example: the 3rd user, female, 168cm, 28 years old can measure as follow:

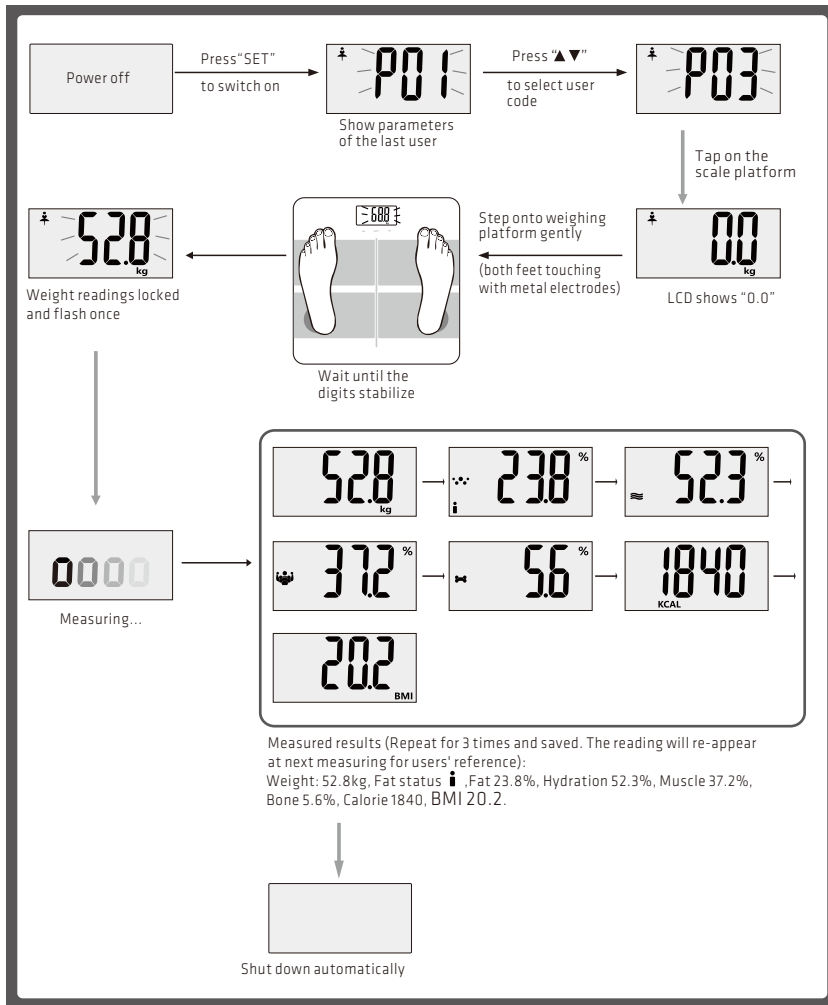TO RECALL TEST RECORD

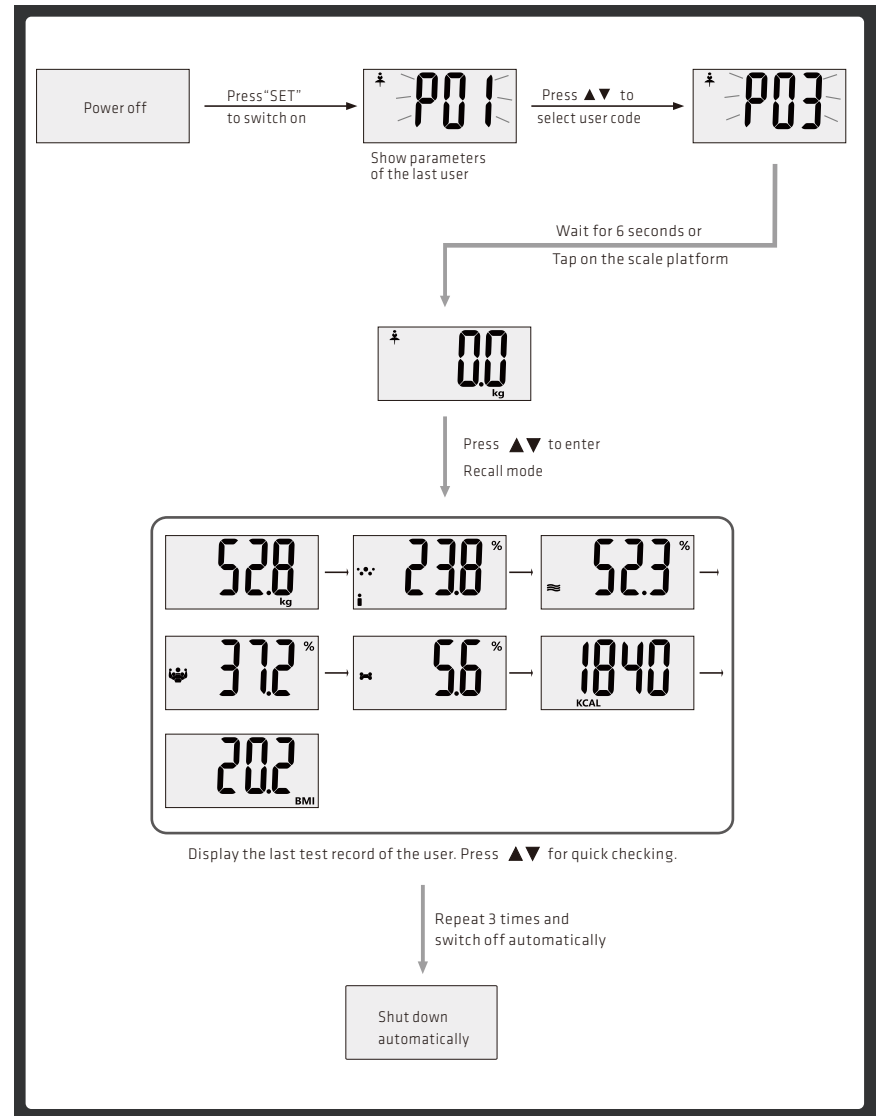WARNING INDICATIONS

Low Battery Indication

The battery power is running low, please replace with a new battery.

Over-load Indication

The weight on the platform exceeds the scale's maximum capacity. Please step off immediately to avoid damage.

Low Fat% Indication

The Fat% is too low. Please lead a more nutritional diet and take good care of yourself.

High Fat% Indication

The Fat% is too high. Please watch your diet and do more exercise.

Re-test

Error exists, please re-test to get the correct result.

ABOUT FAT

1. About Fat%

- Fat% is an index % of fat content in human body.
- The fat content in the human body is of a certain level, overabundance of fat will bring bad effects on metabolism, thereby, various illnesses may occur.

2. Working Principles

As electrical signals of a certain frequency transmit through the human body, the impedance of fat is stronger than that of muscle or other human tissues. Through the use of biological impedance analysis, this product will pass a safe frequency of low voltage electrical current through human body, thus, the Fat% will be reflected according to the resistance differences when the electrical signals are transmitting through human body.

3. Fat/Hydration Content Fitness Assessment Chart

After weighing yourself, you can use the following Assessment Chart for reference:

		Fat%				
Gender	Age	Thin	Slim	Normal	Heavy	Overweight
Male	10-17	< 6%	< 9%	< 26%	< 30%	≥30%
	18-39	< 6%	< 11%	< 22%	< 27%	≥27%
	40-59	< 7%	< 12%	< 23%	< 28%	≥28%
	≥60	< 9%	< 14%	< 25%	< 30%	≥30%
Female	10-17	< 10%	< 20%	< 37%	< 41%	≥41%
	18-39	< 11%	< 21%	< 35%	< 40%	≥40%
	40-59	< 12%	< 22%	< 36%	< 41%	≥41%
	≥60	< 14%	< 24%	< 37%	< 42%	≥42%

		Hydration%		
Gender		Low	Normal	High
Male		< 50%	50-65%	> 65%
Female		< 45%	45-60%	> 60%

The above information is only for reference

ABOUT CALORIES

1. Calorie & Energy Consumption

Calorie: Energy Unit

The data reading indicates the energy you need to sustain the daily metabolism and ordinary activities required for your weight, height, age and sex.

2. Calories & Losing Weight

Controlling your calorie intake combined with proper exercise is an effective way to lose weight. The theory is very simple: when the daily energy intake does not meet your body's energy needs, stored carbohydrates and fats will break down to provide the necessary energy. This process results in weight loss.

LCD has no response

First check that the batteries are still live, if not, replace them if necessary.
Repairs may only be carried out by our authorized customer service or dealers.

ADVISE FOR USE & CARE

1. DO NOT use the scale in humid, excessively hot, or excessively cold environments.
(Temperature range from 5oC to 35oC)
2. DO NOT use the scale on soft mat, lawn, slippery floor or uneven floor.
3. DO NOT drop, shock or strike the scale as it is a precise instrument.
4. Keep the scale clean. DO NOT use acid cleaning agents.
5. DO NOT use the scale in a location with intense vibrations. Stand firmly on the scale when in use.
6. To ensure accuracy, please test at the same time each day.

SPECIFICATIONS

1. Model: POPS3000
2. 74 x 35.6 mm LCD
3. Low battery / over load indication
4. Capacity: 150 kg / 330 lb
5. Graduation: 100g / 0.2 lb
6. Power: 1x3V lithium cell(CR2032)
7. Uses a high precision strain gauge sensor
8. Multifunction: tests body fat %, body hydration %, body muscle%, body bone% and suggests calorie intake requirements and BMI.
9. Body fat % range: 4.0%~75.0%, measurement division: 0.1%
10. Built-in memory for 12 different users
11. Optional normal weighing mode (bypassing parameter setting): Auto step on weighing and consecutive add-on weighings
12. Low battery and overload indication

This symbol on the product or in the instructions means that your electrical and electronic equipment should be disposed at the end of its life separately from your household waste. There are separate collection systems for recycling in the EU.

Notes for battery disposal

The batteries used with this product contain chemicals that are harmful to the environment. To preserve our environment, dispose of used batteries according to your local laws or regulations. Do not dispose of batteries with normal household waste.

For more information, please contact the local authority or your retailer where you purchased the product.

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