NOOA

PRO PERFORMANCE

Physiq Smart Body Composition Scale

PHYSIQ
User Manual

GB

Thank you for purchasing your new NOOA Physiq smart body composition scale. These operating instructions will help you use it properly and safely. We recommend that you spend some time reading this instruction manual to fully understand all the operational features it offers. Read all the safety instructions carefully before use and keep this instruction manual for future reference.

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SE

SAFETY WARNING


1. Read all the instructions carefully before using this scale and keep them for future reference.
2. This product has delicate electronic components. Therefore interference may occur when it is near devices that transmit radio signals. If errors occur in the display, move the scale away from these signals or remove the batteries briefly and then insert them again.
3. Keep batteries away from children. Do not leave batteries lying around. They could be swallowed by children or pets. In the event of swallowing a battery, contact a doctor immediately.
4. This scale is for household use only. Do not use it for commercial purposes.
5. The glass platform will be slippery when it is wet. Keep it dry!
6. Stand still while being weighed.
7. Do not strike, shock or drop the scale.
8. Treat the scale with care as it is a precision instrument.
9. Clean the scale with a damp cloth and avoid water getting inside the scale. Do not use chemical/abrasive cleaning agents.
10. Keep the scale in a cool place with dry air.
11. Always keep your scale in a horizontal position.
12. If the scale fails to turn on, check to see if the batteries are installed, or if their power is low. Insert or replace with new batteries.

13. If there is an error on the display or you cannot turn off the scale for a long time, take out the batteries for about 3 seconds, then install them again to remove soft-fault. If you cannot resolve the problem, contact your dealer for consultation.
14. It's prohibited for people with pacemakers or other implanted medical devices to use this scale.
15. To ensure the accuracy of the data, it is recommended that each measurement be carried out at the same time of day. Don't use in strong electromagnetic conditions or areas with intense vibration.
16. Stand on the scale steadily and touch the electrodes on the platform with your bare feet. Stepping on the edge of the scale may cause it to tip over and lead to a fall. Keep standing firmly on scale when in use.
17. When weighing after not using the scale for a long time or after moving it, please step on the scale surface first, and wait for it to return to zero. Let it switch off automatically before switching on and weighing.
18. To ensure the accuracy of the body composition measurement, please take off your shoes and socks, remove any items that might affect your weight, wash and dry your hands and feet, and then step on the scale with your bare feet and hold the handle with both hands to take the measurement.
19. Measurement results may be misleading after intensive exercise, excessive dieting, or extreme dehydration.


PRODUCT OVERVIEW

PRODUCT OVERVIEW - FRONT SIDE

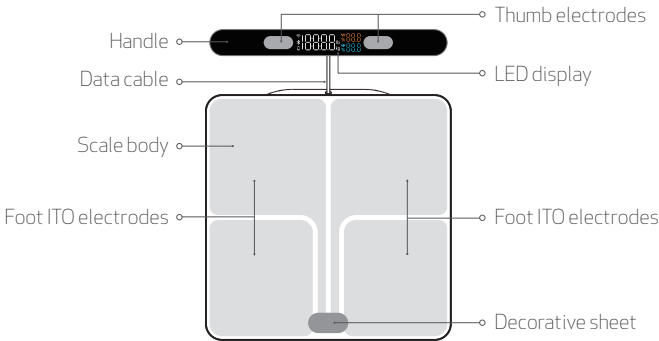PRODUCT OVERVIEW - BACK SIDE

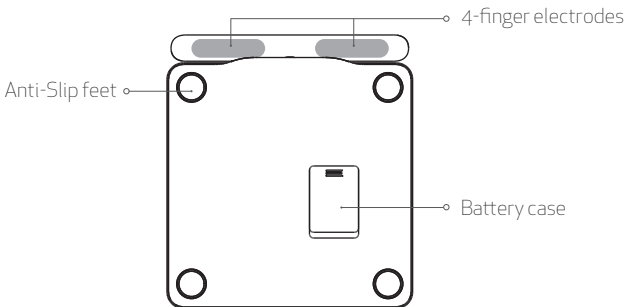BRIEF INTRODUCTION

This scale uses 8 electrodes, multi-frequency BIA and other technologies, which can accurately measure and analyse a full range of body composition data. Used in conjunction with the APP, it can display and record comprehensive body data, as well as analyse and track changes over time. It's the first step to monitoring your weight and maintaining your desired body shape.

IMPORTANT INSTRUCTIONS

The measured data may have some deviations for some individuals, such as:

- Minors, seniors over 80 years old.
- Professional athletes.
- Pregnant women, artificial dialysis patients, people with edema.


MODE INSTRUCTIONS

- Weight measurement mode: No need to pull up the handle. Step on the scale to complete the measurement.
- Accurate body fat measurement mode: You need to pull up the handle according to the correct operation requirements. (For the first time, you need to bind the APP and establish a user profile.)
- Offline Mode: After the first binding, no need to connect the mobile phone. You can save each measurement data, and it will be automatically uploaded and saved to the APP when your mobile phone is connected next time.
- Visitor mode: create a visitor mode through APP, and fill in the information accordingly. Please note that the data in visitor mode will not be saved.
- Small object weighing mode: The scale can directly weigh objects over 0.1kg. For small objects under 10kg, you need to step on the scale surface first to wake up the scale (apply a force greater than 5kg and press the scale surface) and wait until the scale has returned to zero before weighing.

BATTERY INSTALLATION

4 x 1.5V AA batteries are used for this scale. Please replace them as per the instructions below:

Change battery: Open the battery cover, take out the used batteries, then install the new ones properly as shown.

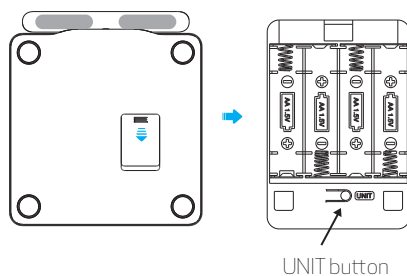WEIGHT UNIT SETTING

Please set the weight unit by APP. If unavailable, you can press UNIT button on the back of the scale to convert the weight unit.

Note:


Please use the correct battery type. Remove batteries if the scale will be unused for more than 3 months.

Usage Tip:

After installing new batteries or moving the scale, it may show a residual value (showing weight when there is none). Please wait about 6 seconds for the scale to switch off and then step lightly on the scale to wake it up. Use it only after it displays 0.00kg again.

OPERATION

Step 1: Place the scale properly and download the APP.

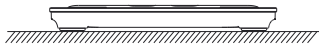Place the scale on hard and flat floor.


Search "MovingLife" in "App Store" or "Google Play" to download the application. Or, for a more convenient method, scan the QR code adjacent.

Step 2: Register and bind the device

Turn on the Bluetooth on your phone, open the APP, and follow the guidelines to register and bind the device:

The APP will prompt you to connect with Wi-Fi after binding. Please follow the APP guidance accordingly.

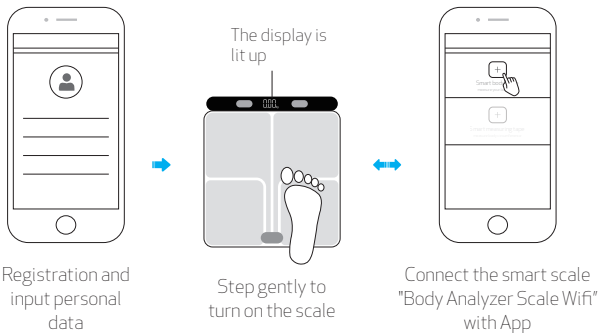Reminder:

Personal body information is used solely for body composition analysis, not for other purposes. To ensure the accuracy of the results, please enter accurate personal information.


Step 3: Measure

• Body Composition Analysis Measurements


1. Place the scale on a hard, flat surface and open the App.
2. Gently step on the scale with bare feet and place your feet on the 4 electrodes on the left and right side of the scale, then lift up the handle. Please hold the handle with both hands as shown in the picture, adjust the height of both hands and straighten your arms (it is best if your arms don't touch your waist).


Step on the scale with bare feet and place your feet on the 4 electrodes

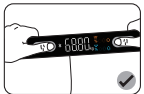
Hold the two front electrodes with thumbs

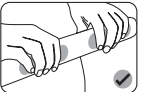
Hold the two back electrodes with the rest of the fingers

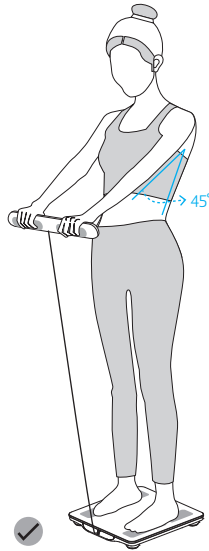• Wrong Measuring Postures

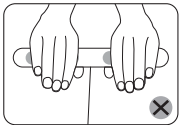
The ten fingers are away from the electrodes

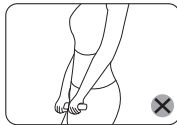
The two arms are not straight or stretched


Squatting down to measure

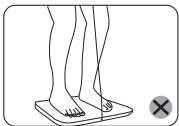The knees are bent or not straightened

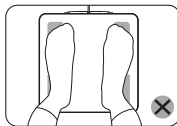Measuring with socks

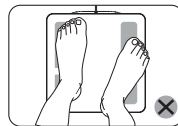The two feet fail to step on the 4 electrodes evenly

3. After the weight value displayed in the window has stabilised and is blinking, the body composition analysis measurement begins.

Body weight measurement

Body fat, body muscle, etc. measuring

Measurement complete

4. Check the measurement reports in the APP

View the measurement reports in the APP and save the historical data.

Tips:

- The whole process will take about 20 seconds. If the APP is not connected when using the scale for the first time, only weight can be measured, not body fat. Please connect the APP before using.
- The APP is regularly updated, so it may differ from what is shown in the manual. For the most current and accurate information, please check the App.

DISPLAY ICONS

	Fat		Muscle
	Wifi		Bluetooth
	Low Battery		%
	Weight Unit		Weight Unit

	Measuring
	Wrong Measurement Indication
	Low Battery Indication
	Over-load Indication
	Too low body fat rate
	Too high body fat rate

MULTIPLE USERS SUPPORTED

App supports registering endless accounts with different phone number or email addresses. Each account can add 12 users at maximum. A maximum of 15 groups of historical data for each user can be stored. Data will be automatically synchronized to the smart device when the APP is running.


Note:

Please make sure only one device is setting up connection with the scale at a time.

If the second device needs to bind the Smart scale, the first device needs to unbind the Smart scale in "Me-My device".

VISITOR MODE

Enter visitor mode in the interface as shown, and fill in the information accordingly, then you can start measuring. Please note that the measuring data in visitor mode will not be saved. If you have already logged in the App with an account, please log out first and press the switch button to switch between visitor mode and normal user mode.

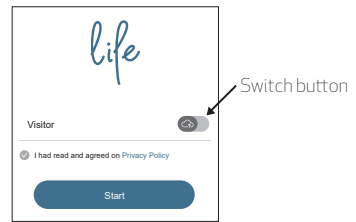FAQ

The software can't obtain weighing results from scale.

- Please confirm the connection icon on the smart device to know the connection status, and check if the WiFi devices are switched on.
- Confirm the Bluetooth is on in your device and the scale is on within the valid range of signal coverage.
- Check if the scale's battery is low; replace if necessary.
- Relaunch the application.
- Turn Bluetooth off and on again, then re-establish the Bluetooth connection between scale and smart phone.

Unable to connect to your device

- Ensure only one device is connecting to the scale at a time. To connect a second device, the first device needs to unbind from the Smart scale in "Me-My device".
- When using the scale with the App open, it will connect via Bluetooth. Without the App open, the scale will use WiFi.
- Please make sure your device supports Bluetooth 4.0 or above.

The screen is unresponsive

- Please first check the batteries and replace them if necessary.
- Repairs may only be carried out by our authorized customer service or dealers.

SPECIFICATIONS

1. Model: PHYSIQ
2. LED digits height 18 mm
3. Low battery / over load indication
4. Capacity: 180 kg / 400 lb
5. Graduation: $100g \leq m < 10kg$ $d = 10g / 0.02lb$
 $10kg \leq m \leq 180kg$ $d = 50g / 0.1lb$
6. Power: 4 x 1.5V AA batteries (not included)

This symbol on the product or in the instructions means that your electrical and electronic equipment should be disposed at the end of its life separately from your household waste. There are separate collection systems for recycling in the EU.

Notes for battery disposal

The batteries used with this product contain chemicals that are harmful to the environment. To preserve our environment, dispose of used batteries according to your local laws or regulations. Do not dispose of batteries with normal household waste.

For more information, please contact the local authority or your retailer where you purchased the product.

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