The logo for 'point' is displayed in white lowercase letters on a red rectangular background with rounded corners.

ASPER NECK & SHOULDER MASSAGER

ASPER
User Manual

The logo 'LIFE MADE EASIER' is enclosed in a red rectangular border. 'LIFE' is in large grey letters, 'MADE' is in smaller orange letters, and 'EASIER' is in large red letters.

GB

Thank you for purchasing your new ASPER NECK & SHOULDER MASSAGER. These operating instructions will help you use it properly and safely. We recommend that you spend some time reading this instruction manual to fully understand all the operational features it offers. Read all the safety instructions carefully before use and keep this instruction manual for future reference.

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SE

IMPORTANT INSTRUCTIONS



READ AND SAVE THESE IMPORTANT SAFETY INSTRUCTIONS

1. **Warning:** This appliance is NOT A TOY. Keep away from infants, babies and children.
2. This appliance can be used by children aged 8 years and older and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they are supervised or have been instructed on the safe use of the appliance and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be carried out by children without supervision.
3. DO NOT use this appliance for more than 20 minutes.
4. DO NOT fall asleep while using this appliance.
5. Pregnant women, infants, children, people with a swollen or inflamed body parts, people with serious muscle or joint aches, and anyone wearing a pacemaker, are not recommend to use this appliance.
6. DO NOT operate this appliance after taking pain killers, sedatives or alcohol.
7. Close supervision is required when disabled or ill individuals operate this appliance.
8. Massage should be comfortable and pleasant, if you feel pain or discomfort, discontinue use and consult your doctor.
9. DO NOT place the appliance near water, under direct sunlight or in excessive heat.



10. Do not use when lying flat or reclining.
11. Do not use this product during precise tasks or driving.
12. Store out of children's reach.
13. Children under 8 should not use this product. Those aged 8-12 should only use it under adult supervision.
14. Do not scratch or puncture this product with sharp objects.
15. Do not disassemble the internal components of the product to prevent damage.
16. This is a household electronic healthcare product, not a medical device, and not suitable for diagnosis or treatment.
17. Remove jewelry from the area being massaged before use.
18. Use the provided charging cable and a standard-compliant adapter, ensuring safe electrical use.
19. Excessive use may cause discomfort or skin irritation.

Product Overview

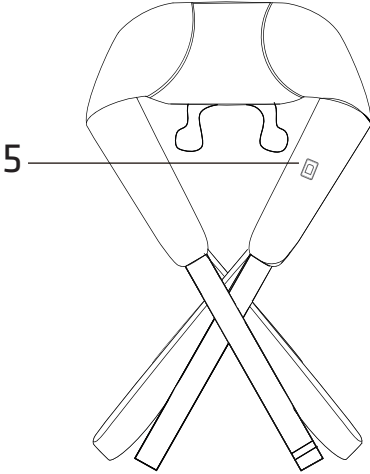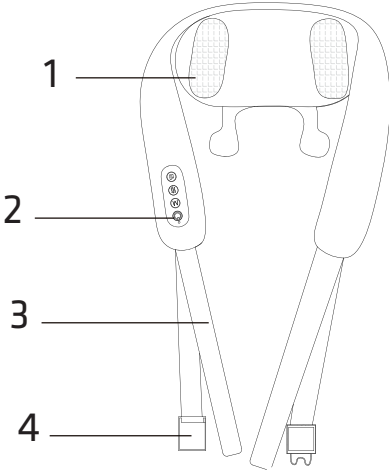- 1. Massage node
- 2. Control panel
- 3. Hand strap

- 4. Elastic strap buckle
- 5. TYPE-C charging port

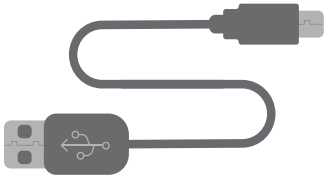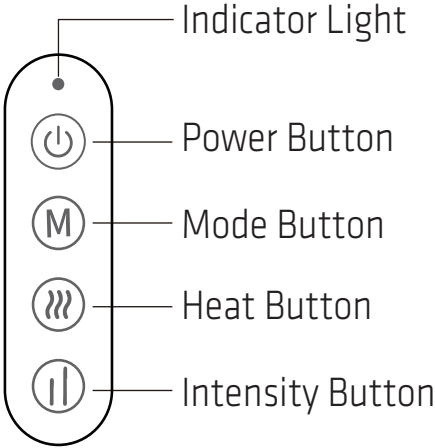USB-C Charging Cable

Operating Instructions

Power button

Power on/off

1.Power on: Press and hold this button for 1 seconds, the default mode of power-on is Mode 1, the hot compress is turned on, the intensity is Gear 1.

2.Power off: Press and hold for 1 seconds again to turn off the massager.

Mode key

Switch mode

Short press to switch the massage mode in a cycle.There are 2 massage modes in total, namely: mode 1 (vitality mode) and mode 2 (relaxation mode).

Heating key

Hot compress mode


Short press the button to activate or deactivate the hot compress function. The indicator light will display solid red when the function is active, and switch to green when deactivated.

Intensity key

Massage intensity switching

Short press the button to cycle through two massage intensity levels.


Massage Position

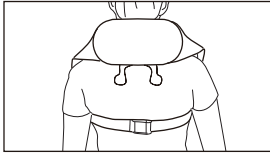Hands-free buckle strap design

Neck massage

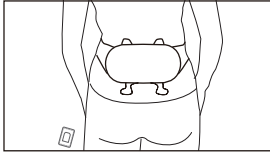
Lower back relaxation

Relief for leg pressure

Preparation before use

- Visually examine the product for shipping damage prior to initial operation.
- Charge the product for about 2.5 hours to fully charge it before the first use.
- Connect the USB-C cable to charging port of the product, then to an certified adapter(5V 2A), plug the adapter to an AC power outlet.

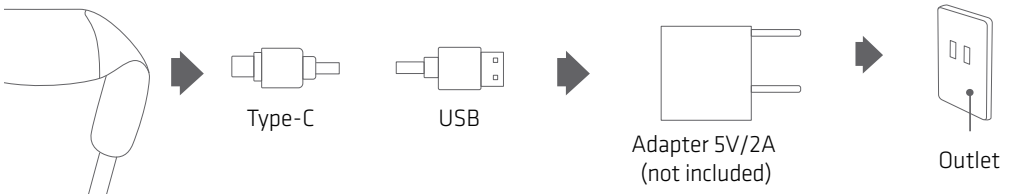- While charging, the product cannot be turned on and used.
- The battery indicator light is red during charging and turns green when fully charged.
- If the battery level drops below 20%, the indicator flashes red to signal low battery; charge it promptly

- Remove necklaces, metal accessories, etc. before use to avoid discomfort.
- Adopt a comfortable posture and completely relax your body.
- Wear the massager on your shoulders and neck, and adjust it to the most comfortable position, and then start the massage.

Note: It is normal for the device and the power adapter to heat up slightly during charging.

Care and Maintenance

1. Keep the product clean and away from places with high humidity or high temperature and environments with volatile gases.
2. Protect against moisture, rodents, and dirt; do not place heavy items on the product.
3. In case of failure, seek professional repair services. Non-professionals should not attempt disassembly.
4. Stop using the product if you notice any of the following: a burnt smell, intermittent power disconnections, a hot power cord, cord damage, or any other abnormalities.

Troubleshooting

This section summarizes the most common problems you could encounter with the product. If you are unable to solve the problem with the information below, contact the service center immediately.

Problem	Possible cause	Solution
Can not power on	Battery power is insufficient.	Charge it, then try again.
The product does not work	Low battery	Charge it, then try again.
	Extended use, the motor's overheating protection has activated.	After the motor has cooled down naturally, restart it and continue to use.
Stops during use	The 10-minute default working time is up.	Long press the power button to restart.
	Low battery	Charge it and then use it again
Weak massage	Low battery	Fully charge it.
No heat compress	The hot compress function is not turned on	Short press the heating button to turn on the hot compress function.
Internal odor or smoke during operation	Product failure	Stop using and turn of the power supply, then contact the service center.

Problem	Possible cause	Solution
The charging indicator light is not on	The USB-C socket has poor contact /poor power contact.	Pull and plug the Type-C charging cable / power cord.

Specification:

Model: ASPER

Input: 5V \Rightarrow 2A

Rated power: 10W

GB

This symbol on the product or in the instructions means that your electrical and electronic equipment should be disposed at the end of its life separately from your household waste. There are separate collection systems for recycling in the EU. For more information, please contact the local authority or your retailer where you purchase the product.

Power International AS
PO Box 523, N-1471 Lørenskog, Norway

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