point

FELIX NECK & SHOULDER MASSAGER

FELIX
User Manual

LIFE MADE EASIER

GB

Thank you for purchasing your new POINT FELIX NECK & SHOULDER MASSAGER. These operating instructions will help you use it properly and safely. We recommend that you spend some time reading this instruction manual to fully understand all the operational features it offers. Read all the safety instructions carefully before use and keep this instruction manual for future reference.

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IMPORTANT SAFETY INSTRUCTION

When using the product, all instructions and basic precautions should always be followed, including the following:

DANGER - To reduce the risk of electric shock: Always unplug this product from the power outlet immediately after charging and before cleaning.

WARNING - To reduce the risk of burns, fire, electric shock, or injury to persons:

1. The product should never be left unattended when plugged in.
2. Do not operate under blanket or pillow to prevent excessive heating that can cause fire, electric shock, or injury to persons.
3. Ensure air openings are not blocked during operation.
4. Keep the air openings free of lint, hair, and similar debris.

5. Do not drop or insert object into any openings.
6. Do not use outdoors.
7. Do not operate where aerosol (spray) products are being used or where oxygen is administered.
8. Do not plug or unplug the product with wet hands.
9. Do not carry this product by its supply cord or use the cord as a handle.
10. Remove the power plug from the socket immediately after a power failure to prevent electric shock, short circuits, or fire.
11. Never operate this appliance if it has a damaged cord or plug, is not working properly, has been dropped or damaged, or has fallen into water. Return it to a service center.
12. Use this appliance only for its

intended use as described in this manual.

13. Keep the cord away from heated surfaces. If any fault is found, stop using immediately and cut off the power supply to avoid accidents or injuries. Do not attempt to repair it by yourself; contact the service center immediately.
14. For charging, use only a protection class II power supply approved for household products, with SELV output 5V, Max. 2A.
15. Charge the appliance only at the safety extra low voltage specified on the appliance.
16. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and

knowledge, unless they are supervised or instructed on the use of the appliance by a person responsible for their safety.

17. Children should be supervised to ensure that they do not play with the appliance.
18. This appliance contains non-replaceable batteries.
19. The appliance has a heated surface; users insensitive to heat should exercise caution.

Check that the product and accessories are intact and undamaged before use.

IF YOU FEEL ANY ABNORMAL SYMPTOMS DURING USE, STOP IMMEDIATELY AND CONSULT YOUR DOCTOR. CONTINUED USE MAY BE HARMFUL TO YOUR HEALTH.

This product is not recommended for the following individuals.

1. Those susceptible to electromagnetic malfunctions, such as individuals with cardiac pacemakers and implantable medical devices.
2. Individuals with impaired temperature sensation.
3. People with malignant tumors.
4. People with heart disease.
5. Individuals who are pregnant or have recently given birth.

6. Individuals with an abnormal or curved spine.
7. Those with sensory disturbances due to peripheral circulation issues in diabetes.
8. People with scars on areas where the massager will be applied.
9. Individuals who require quiet environments or are in poor health.
10. Those with osteoporosis, severe fractures, sprains, muscle pain, or acute head pain.
11. Individuals with a fever (38°C or higher).
12. Those experiencing acute inflammatory symptoms, fatigue, chills, blood pressure changes, or intense stress.

13. People with tendonitis.

14. Individuals receiving treatment at medical facilities.

The device is not to be used for commercial use or medical purposes. If you have health concerns, consult your doctor before using this product.

Only use indoors!


Do not use the device in damp rooms (e.g. when bathing or showering)


- Do not wash!
- Do not dry clean!
- After using this product for 10 minutes, rest for 30 minutes before use again.

PRODUCT OVERVIEW

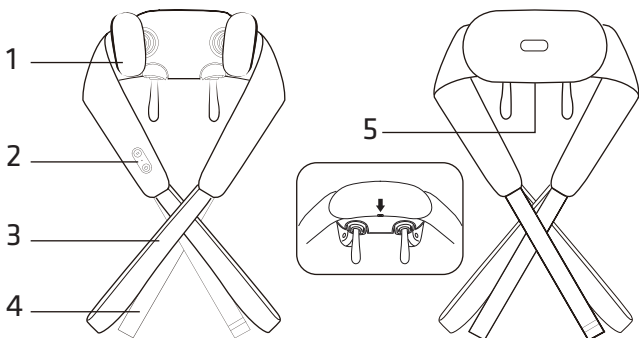
1. Massage node
2. Control panel
3. Hand strap

4. Elastic strap buckle
5. TYPE-C charging port
(on the bottom)

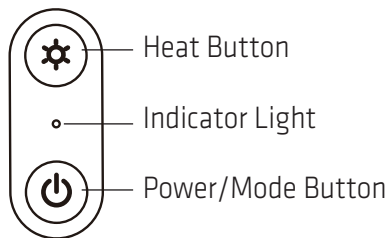USB-C
Charging Cable

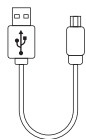GB


This product is suitable for home use only

* Use only certified chargers when charging this device.

This product comes with a USB-C 5V charging cable. To charge the product, you need to use a compatible adapter. It is recommended to use a 5V/2A adapter for optimal performance and safety.

* If the product is not used for an extended period, it is recommended to charge it once every 6 months to maintain battery health and prevent degradation.


Massage Position

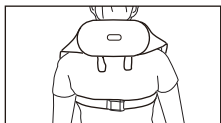Hands-free buckle strap design

Neck massage

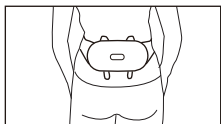Lower back relaxation

Relieving leg pressure

OPERATING INSTRUCTIONS

1. Power On: Long press the power button for 2 seconds to start the machine. The default mode is Kneading, indicated by a white light.
2. Mode Switching: Press the mode button to cycle through the two modes: Kneading → Shiatsu → Kneading.
3. Power Off: Long press the power button for 2 seconds to turn off the machine.
4. Automatic Shutdown: The machine will automatically shut down after 10 minutes of continuous use.
5. Heat Function
 - By default, the machine starts without heat.
 - Press the Heat Button to activate heat; the light turns red. Press again to turn off heat; the light returns to white.

HOW TO CHARGE

1. The machine cannot operate while charging.
2. When the machine is charging, the indicator light turns red.

- When the battery is fully charged, the indicator light turns white.
- If battery power drops below 20%, the indicator flashes red, signaling it's time to charge promptly.


Note: It is normal for the unit and charging plug to slightly heat up during charging.

Charging Time

The charging process takes approximately 3 hours.

After Charging

Once charging is complete, unplug the power plug from the outlet.

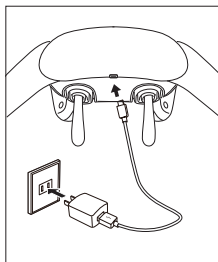MAINTENANCE

Cleaning


Before cleaning, ensure the power is turned off and the power plug is disconnected from the socket.

- Remove dirt by wiping with a cloth moistened with soapy water. Ensure the cloth is thoroughly wrung out prior to use.

- Avoid the use of alcohol, paint thinner, benzine, or nail polish remover, as these substances may cause damage to the appliance, including cracking or discoloration of its components.

Storage

Keep the product clean and far from the wet environment.

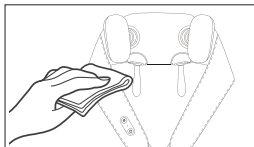GB

Troubleshooting

| Problem | Cause | Remedy |
|----------------|---|---|
| Can not start. | The power button is not on. | Long press the power button to turn on. |
| | 10-minute automatic massage time is up. | Long press the power button to turn the massager on again . |
| | An excessive load applied to the product. | Reduce the load and start again. |

| Problem | Cause | Remedy |
|--------------------------------------|---|--|
| Sudden failure to start during use. | The product's battery is fully drained. | Charge promptly. |
| Noise can be heard from the product. | Noise is made by the mechanical parts of the product. | No action needed. |
| The product is smoking. | Mechanical issue or damage. | Stop using immediately and contact the service center. |

If you do not use it for a long time, ensure that the function is normal before using it again.

Periodical inspection.

Check for any of the following faults:

- Power cord is damaged, deformed, or smells burnt.
- Abnormal sounds or vibrations during operation.
- Product does not work properly after being switched on.
- Leather case or appearance is damaged.

If any of these faults occur, stop using the product immediately and disconnect from the power supply. Send it to a designated maintenance point for repair.

GB

If you still experience problems after following these methods, please contact our customer support team.

SPECIFICATIONS

Model: FELIX

INPUT: DC 5V $\overline{=}$ 2A

Power consumption: 18 W

Built-in Battery: 7.4V/2000mAh lithium battery

This symbol on the product or in the instructions means that your electrical and electronic equipment should be disposed at the end of its life separately from your household waste. There are separate collection systems for recycling in the EU. For more information, please contact the local authority or your retailer where you purchased the product.

Power International AS
PO Box 523, N-1471 Lørenskog, Norway

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