**point**

# CURE BACK MASSAGER

CURE  
User Manual

**LIFE MADE EASIER**

## GB

Thank you for purchasing your new POINT CURE BACK MASSAGER. These operating instructions will help you use it properly and safely. We recommend that you spend some time reading this instruction manual to fully understand all the operational features it offers. Read all the safety instructions carefully before use and keep this instruction manual for future reference.

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
## **SAFETY INSTRUCTIONS**

1. People with medical conditions such as pacemakers, heart disease, tumors, spinal issues, osteoporosis, skin infections, recent surgeries, or those who are pregnant, post-delivery, or undergoing medical treatment should consult a physician before using this product.
2. Do not use the product in wet environments like bathrooms or near pools. Ensure no water spills on or enters the product.
3. Do not use the product in environments with high humidity or temperature or near heat sources like heaters.
4. Charge only at the specified low voltage indicated on the appliance.
5. To prevent electrical hazards, do not open the product's outer cover.
6. This product is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, unless they are supervised or instructed on the use of the product by a person responsible for their safety. Children should be supervised to ensure that they do not play with the product.
7. Always return the product to a service center for examination or repair if the product is broken; do not attempt to repair the product yourself.




8. Only position the product as described in this manual.
9. Do not use the product on a chair with faulty or unstable back support.
10. Do not apply medical or cosmetic creams or lotions on the product or on your back before or during the use of the product.
11. Do not shake or drop the product.
12. Ensure the product is not placed under heavy objects or squeezed between items during movement to prevent damage.

# PRODUCT OVERVIEW

Front

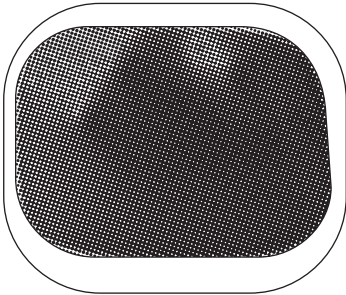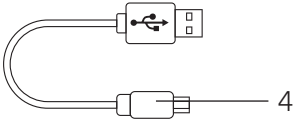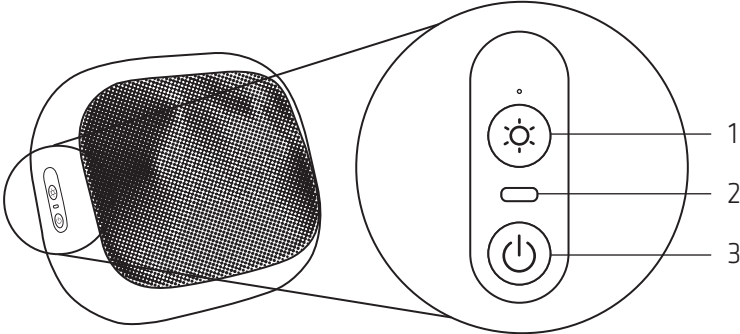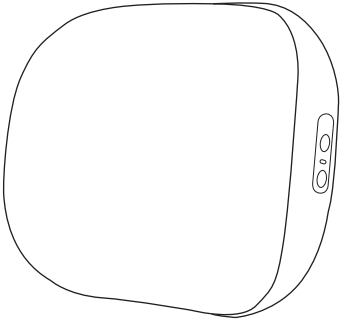Back

- 1.Heat & vibration switch for back side
- 2.USB charging port
- 3.Massage switch for front side
- 4.USB-C charging cable

# OPERATING INSTRUCTIONS

## How to operate

1. Press  and hold for 2 seconds. The product will start its kneading massage and heat function for 10 minutes by default.
2. Short press  to switch the kneading speed between level 1 (green light) and level 2 (red light).
3. Press  and hold for 2 seconds to turn off the product.

## Heating & vibration

1. When the product is power off, press  and hold for 2 seconds to turn on the heat/ vibration function on the back side; by default, it starts with just heat (red light).
2. Short press to switch from vibration only (green light), heat + vibration (orange light), heat only (red light).

**Note:** The default usage time for this product is 10 minutes. Do not use it for more than 30 minutes at a time.

## Charging

1. Do not use the product while charging.
2. If power drops below 20%, the red light will flash quickly. Connect the device to the charger (DC 5V 2A) immediately.
3. During charging, the red light will flash slowly.
4. Once fully charged, the light will turn green and stay on constantly. Disconnect from the charger.


## Application:


Alleviates stiff neck

Soothes back fatigue

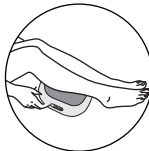Relieves leg fatigue

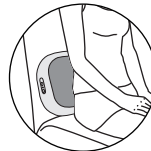Relieves lower back

Eases menstrual discomfort (back side)

## TROUBLESHOOTING

This section summarizes the most common problems you could encounter with the product. If you are unable to solve the problem with the information below, contact the service center immediately.

<b>Problem</b>	<b>Possible Cause</b>	<b>Solution</b>
The product does not work	The product has low battery.	Recharge the product until the battery is full.
	The product has just finished its 10-minute massage and has automatically switched to standby mode.	If you want another massage, long press the Power Switch to turn it on again.
	You may have over-used the product, causing the overheat protection mechanism to switch off the product.	Unplug the product and let it cool down for 30 minutes before you use it again.
	An excessive weight has been put on the product (e.g. someone has stood on the product or a heavy object has been put on top of it).	Remove the excessive weight on the product.
The product is warm	You may have overused the product / turned on the heat unintentionally.	Let it cool for 30 minutes/ Turn off the heat function.

## CARE & MAINTENANCE

### Cleaning

- 1.Never use scouring pads, abrasive cleaning agents or aggressive liquids such as petrol or acetone to clean the product.
- 2.Do not immerse the product in water.
- 3.Clean the product with a dry cloth.

### Storage

- 1.Store the product in a cool and dry place, and keep it away from sunlight and humidity.
- 2.Do not put heavy objects on top of the product.

## SPECIFICATIONS

Model: CURE

Input: DC 5V  $\overline{=}$  2A

Battery: 2200mAh lithium battery

Charging time: Approx. 3.5hr

GB

This symbol on the product or in the instructions means that your electrical and electronic equipment should be disposed at the end of its life separately from your household waste. There are separate collection systems for recycling in the EU. For more information, please contact the local authority or your retailer where you purchase the product.

Power International AS  
PO Box 523, N-1471 Lørenskog, Norway

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