NOOA

PRO PERFORMANCE

# Body Analyser Scale

NO902  
User Manual

GB

Thank you for purchasing your new NOOA body analyser scale. These operating instructions will help you use it properly and safely. We recommend that you spend some time reading this instruction manual in order that you fully understand all the operational features it offers. Read all the safety instructions carefully before use and keep this instruction manual for future reference.

SE

Tack för att du köpt en ny NOOA kroppsanalysvåg. These operating instructions will help you use it properly and safely. We recommend that you spend some time reading this instruction manual in order that you fully understand all the operational features it offers. Read all the safety instructions carefully before use and keep this instruction manual for future reference.

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## SAFETY WARNING

1. Read all the instructions carefully before using this attachment and keep them for future reference.
2. This device has delicate electronic components. Therefore interference may occur when it is near devices which transmit radio signals. If errors occur in the display, remove away from these signals or remove the battery briefly and then insert it again.
3. Keep batteries away from children. Do not leave batteries lying around. They could be swallowed by children or pets. In the event of swallowing a battery, contact a doctor immediately.
4. This scale is for household use only. Do not use for commercial purpose.
5. The glass platform will be slippery when it is wet. Keep it dry!
6. Keep standing still during the course of weighing.
7. Do not strike, shock or drop the scale.
8. Treat the scale with care as it is a precise instrument.
9. Clean the scale with a damp cloth and avoid water getting inside the scale. Do not use chemical/abrasive cleaning agents.

10. Keep the scale in cool place with dry air.
11. Always keep your scale in horizontal position.
12. If failed to turn on the scale, check if battery is installed, or battery power is low. Insert or replace new battery.
13. If there is error on display or can not turn off the scale for a long time, please take off the battery for about 3 seconds, then install it again to remove soft-fault. If you can not resolve the problem, please contact your dealer for consultation.

## PRODUCT OVERVIEW



## REPLACING BATTERIES

Depending on the model of your scale, the scale may operated by one CR2032 (3V) lithium battery or two CR2032 (3V) lithium batteries (included). Before first use, remove the insulation strip at the battery compartment.

1. Open the battery compartment cover.
2. Remove the old batteries (if any) using a sharp object as shown in the picture. Install one CR2032 (3V) lithium battery or two CR2032 (3V) lithium batteries. Make sure the positive (+) is facing upwards.

HINT: Make sure one side of the battery is placed underneath the battery baffle, then push the other side of the battery until it clicks and secured.

3. Close the battery compartment cover.

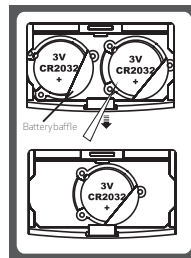
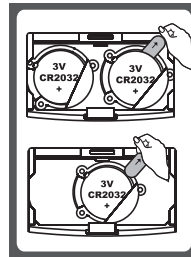
## PERFORMANCE INDEX

Max Weighing Capacity	Min Weighing Capacity	Division	Fat% Range	Fat% Division	Hydration% Range	Hydration% Division	Muscle% Range	Muscle Division
150kg 330lb	2kg 4.4lb	0.1kg 0.2lb	4.0% -60.0%	0.1%	27.5% -66.0%	0.1%	20.0% -56.0%	0.1%

## CHANGING WEIGHT UNIT

Depending on the model of your scale, your scale may support weight units in kg/lb/st.


Press the UNIT button repeatedly at the rear of the scale to select your desired unit of weight. The units will be displayed in the LC display.

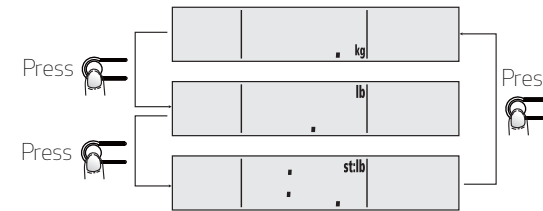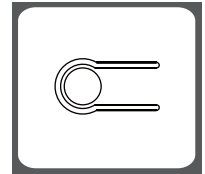After the weight unit is confirmed, the scale will automatically switch off in 3 seconds. The selected unit will be shown on the display when the scale is switched on again.

HINT: If the weight unit "lb/st" is selected, "feet/inch" is the default height unit; if the weight unit "kg" is selected, "cm" is the default height unit.

## LC DISPLAY
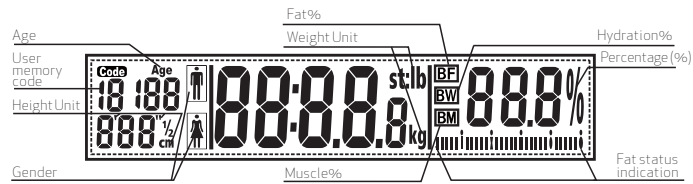
## FUNCTION TOUCH BUTTONS


- ⊙ SET Press to switch the scale on.  
Press to confirm the value entered in setting personal profile.
- ▲ Press the Up/ Down button to select your desired value in setting personal profile.
- ▼ Press the Up/ Down button to select your desired value in setting personal profile.

HINT: Press and hold the Up/ Down button to speed up the process of setting the values.

## NORMAL WEIGHING MODE (without entering personal data)

### (1) Auto step on weighing

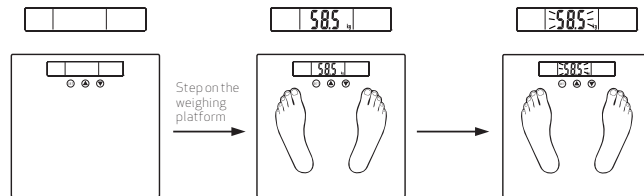
- Place the scales on a solid, flat horizontal surface. Do not place on a carpet or any soft surface. Always make sure there are no other objects beneath the scale before weighing.
- Place even weight on each leg and do not move. Moving during the weighing process can lead to discrepancies in the results. The scales will start to measure immediately. After flashing, the weight display will show the result.
- Step off the platform and view the weight measurement.

### (2) Consecutive add-on weighing

Consecutive add-on weight can also be measured after the first weight shown on the display.

NOTE: The add-on weight should exceed 2kg.

HINT: This function is useful if you want to weigh your baby. See below procedures for details.

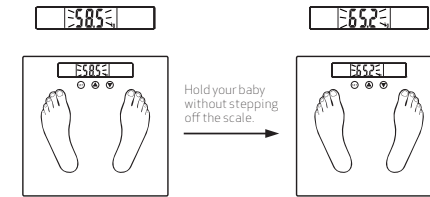- First weigh yourself by placing even weight on each leg and do not move. The scales will start to measure immediately. After flashing, your weight will be shown on the display.
- Do not step off the platform. Now hold your baby and weigh together. The scales will start to measure immediately. After flashing, the new weight measurement will be displayed.
- Step off the platform and view the new weight measurement. By subtracting your own weight from the new weight measurement, you will get the weight of your baby.

## PERSONAL BODY ANALYSIS MODE

### (1) Entering personal data

- First time use, the personal data settings are defaulted as below:

Gender	Height	Age
Male	165cm/5'05.0"	25

- Personal data's range can be set as below:

Memory	Gender	Height	Age
1-12	male/female	100~250cm 3'03.5"-8'02.5"	10~100

## WARNING INDICATIONS


### Low Battery warning

If the scale shows "Lo", it indicates that the batteries are too weak and it is time to replace the batteries. Refer to the section "REPLACING BATTERIES" to replace batteries.


### Overload warning

If the scale shows "Err", it indicates that the scale is overloaded. The maximum capacity is indicated on the scale. Step off the scale to avoid causing damage to the scale.

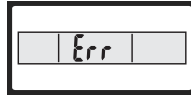
### Low Fat warning

If the scale shows "ErrL", it indicates that your body fat portion is too low. Please consult your doctor or nutritionist for professional advice.

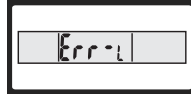
### High Fat warning

If the scale shows "ErrH", it indicates that your body fat portion is too high. Please consult your doctor or nutritionist for professional advice.
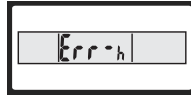
### Error warning

If the scale shows "C", it indicates that there are measuring errors. Please step off the scale and retry.
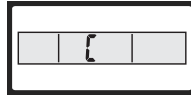
## ABOUT FAT

### 1. About Fat%

Fat % is a % index of fat content inside a human body. The fat content in a human body is in a stable level. Any big fluctuations in fat will have a bad effects on the metabolism, which may even cause various illnesses. Hence it is always very important to keep track on our fat% in order to have a healthy body.

### 2. Working Principles

By applying a low-voltage electrical current of a safe frequency through the human body, with the fact that the impedance of fat is bigger than that of muscle or other body issues, this body analyser scale can analyse the biological resistance which reflect the fat% in the body.

## 3. Fat/Hydration Content Fitness Assessment Chart

After attaining the analysed result, you may use the below assessment table for reference.

Age	Fat(F)	Hydration(F)	Fat(M)	Hydration(M)	Status	Body fat index
≤30	4.0-16.0%	66.0-57.8%	4.0-11.0%	66.0-61.2%	Underweight	
	16.1-20.5%	57.7-54.7%	11.1-15.5%	61.1-58.1%	Thin	
	20.6-25.0%	54.6-51.6%	15.6-20.0%	58.0-55.0%	Normal	
	25.1-30.5%	51.5-47.8%	20.1-24.5%	54.9-51.9%	Fat	
	30.6-60.0%	47.7-27.5%	24.6-60.0%	51.8-27.5%	Overweight	
>30	4.0-20.0%	66.0-55.0%	4.0-15.0%	66.0-58.4%	Underweight	
	20.1-25.0%	54.9-51.6%	15.1-19.5%	58.3-55.3%	Thin	
	25.1-30.0%	51.5-48.1%	19.6-24.0%	55.2-52.3%	Normal	
	30.1-35.0%	48.0-44.7%	24.1-28.5%	52.2-49.2%	Fat	
	35.1-60.0%	44.6-27.5%	28.6-60.0%	49.1-27.5%	Overweight	

- The above information is for reference only.

## ⚠ CARE AND MAINTENANCE

- Clean the scale using a slightly moist cloth. Do not use chemical cleaners or abrasive cleaning agents for cleaning the scale. Ensure there is no water entering the scale. If this does occur, then wait until the device has dried out fully before you use it again.
- Always place the scales on a solid, flat horizontal surface. Do not place on a carpet or any soft surface.

## Features/Specifications

- 1.Size of LCD: 160x30mm
- 2.Using high precision strain gauge sensor
- 3.Touch button
- 4.Auto step on weighing and consecutive add-on weighing without entering personal data
- 5.Analysing percentage of body fat, body water and body muscles
- 6.Built-in memory for 12 different users
- 7.Analyse body fat condition into 5 levels
- 8.Low battery and overload indication
- 9.Memory recall from last user

This symbol on the product or in the instructions means that your electrical and electronic equipment should be disposed at the end of its life separately from your household waste. There are separate collection systems for recycling in the EU.

### Notes for battery disposal

The batteries used with this product contain chemicals that are harmful to the environment. To preserve our environment, dispose of used batteries according to your local laws or regulations. Do not dispose of batteries with normal household waste.

For more information, please contact the local authority or your retailer where you purchased the product.

Imported and Exclusively marketed by:  
Expert AS, PO Box 43, N-1483 Hagan

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